



NEONATAL

SIRIO BABY 200

Portable ICU ventilator for children and newborns



THAI MEDICAL TRADING CO.,LTD
67/17 Moo.1 Mahasawat,Phutthamonthon
Nakhonpathom 73170 Thailand 081-3451047



Simple ventilation advanced Technology

SIRIO BABY 200

The project of the new **Sirio Baby 200** started in 1982 and now, this friendly user equipment has been refined achieving the goal of combining an exceptional flexibility with a very high level technology.

This ventilator has been conceived especially for the emergency when a ventilatory support is needed during the transport. Solid and light at the same time and equipped with a long duration back-up battery which grants the maximum flexibility.

Versatility

The new **Sirio Baby 200** is able to offer a wide range of therapeutic treatments as well as to perform a pressure and volume control ventilation guaranteed either with an invasive or with a non-invasive modality, able to meet all the respiratory needs of newborns and prematures (weight < 3kgs).

Simplicity

The intuitive and versatile interface carefully studied for doctors operating with neonatal emergency situations makes this ventilator the ideal equipment to be used along with the transport incubator and for the transport in ambulance, mobile clinics, helicopter and airplanes.

Performances

The new **Sirio Baby 200** grants excellent performances taking advantage of Siare's high standards technology. Equipped with different ventilation modalities: CPAP – IMV – IPPV – MANUAL – PEEP – Pressure trigger – FiO₂ adjustable from 21% to 100%. Inspiratory and Expiratory time adjustable. Integrated spirometry. Long duration back-up battery and a wide range of dedicated accessories grant an easier mobility of the device.

MRI Version

The new **Sirio Baby 200** is manufactured also in the amagnetic version (item code 09637/RMN) for the use in the Magnetic Resonance unit (MRI) from 1,5 to 3,0 Tesla. The amagnetic version is the ideal choice for the transport of newborns from the resuscitation room to the magnetic resonance unit with a continue ventilation support granting to the newborns a higher comfort and security.



The new Sirio Baby 200 can be used in Magnetic Resonance (MRI) from 1,5 to 3,0 Tesla.



The new Sirio Baby 200 can be used also in helicopter, mobile clinics, ambulances and airplanes.

Acoustic and visual alarms

Simple and intuitive settings

Built-in spirometer



Adjustable trigger

Long duration battery


Patient circuit with expiratory valve and proximal flow sensor



Transport version with professional bag for ventilator and 2 lt cylinder with inox plate for being hooked/unhooked in the ambulance.

Transport version with professional inox support for ventilator and 2 lt cylinder with support for the relevant installation in the ambulance.

DATA FOR ORDERING

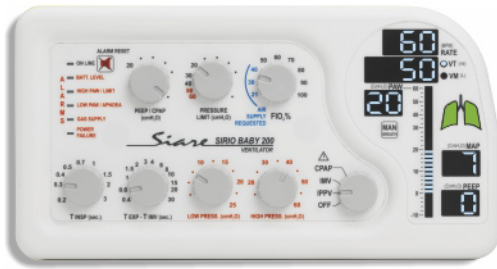









Item Code 09637
SIRIO BABY 200 Neonatal Portable ICU Ventilator

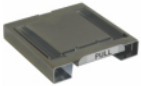

Item Code 09637/RMN
SIRIO BABY 200 Neonatal Portable ICU Ventilator for MRI

OPTIONAL ACCESSORIES

 <p>Code A26.049048 Stand for ventilator with 5 wheels</p>	 <p>Item Code G00203001 Wall plate for rapid connection of stainless steel support</p>	 <p>Item Code 104641/SLN Double-hose reusable Hytel neonatal patient circuit w/exp. valve and flow sensor</p>
 <p>Item Code 05048/RMN Castor LT, Trolley for MRI</p>	 <p>Code G00213000 Plate for rapid connection of ventilator on ambulance flat surfaces</p>	 <p>Item Code A17.07000001 Proximal flow sensor</p>
 <p>Item Code 06750/SB Aluminium supporting arm</p>	 <p>Item Code G90003012 Transport bag for ventilator and cylinder</p>	 <p>Item Code 03000 Vapour 2, Electronic Humidifier</p>
 <p>Item Code 10001/2 2 lt O2 cylinder Item Code 108001/CM O2 pressure reducer</p>	 <p>Item Code G00208003 Tubular support for ventilator and protection strip</p>	 <p>Item Code G80063001 Long duration battery pack</p>

For additional accessories please refer to the price list in force