

# ZOOMAX GOLD VETELEV SYSTEM



## VETELEV ELEVATING TABLE

A SIMPLIFIED AND COST-CONSCIOUS VERSION OF OUR POPULAR PHOENIX MEDICAL PATIENT TABLES, FEATURING



### SUBSTANTIAL FREEDOM OF POSITIONING

- 36 cm vertical movement
- Floating table top with electromagnetic brakes (152 cm x 71 cm with extension option)
- Over 100 kg lifting capacity



### TOTAL SAFETY AND COMFORT

- Integrated crash guard sensors on the table
- Grid cabinet for standard cassette size wifi or fixed DR panels ensuring crystal clear images
- Configuration for fixed or removable grids, including safe „parking slot” for unused grid
- Wireless foot pedal (with long-life battery) to minimize cabling (optional)
- Foot pedal doubles as light field switch



## ZOOMAX TUBE STAND

- Sleek, robust design
- Fixed position tube with 100-136 cm SID (with Vetelev table)
- Free-standing setup for total positioning freedom (optional)



## SELECTION OF CONTROL CONSOLES

- 22" Touchscreen Monitor Console to set exposure parameters and for full image capture and analysis
- 10" LCD console – integrated into X-ray tube cradle – to set exposure parameters, featuring:
  - slim, modern User Interface
  - customizable main page and workflow buttons (e.g. APR, grid control, focus, etc.)
  - remote diagnostics and software updates
  - easy installation
- ZooMax LC or CPI Rad membrane consoles with comfort of all key generator parameter settings, including APRs and advanced controls such as AEC




## INTEGRATED GENERATOR CONTROL AND IMAGE ACQUISITION

Single console seamless workflow:


- Unlimited APRs
- Full DICOM compatibility

Sophisticated image processing features:

- Detector status indicators
- Wifi signal strength indicator
- Wide variety of DR detector choices
- Portable solutions with tablet, laptop workstation
- Mini PACS included
- Cloud storage available
- Mobile viewing on any device

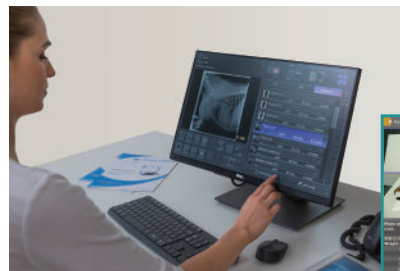

Patient positioning guide

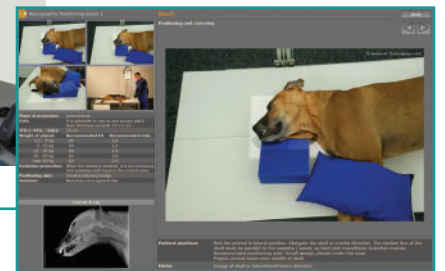

## HIGH QUALITY COMPONENTS


- Industry-leading x-ray tubes
- High-density grids
- Certified manual collimators with light beam or laser for easy positioning
- Powerful generators for minimal exposure times

## GENERATOR OPTIONS

- 32 kW to 50 kW output range
- Stored energy generators up to 50 kW (capacitor discharge) at 15 Amp
- 2.8 kW to 5 kW monoblock generators

## REMOTE SUPPORT

Service software allows troubleshooting to minimize downtime.

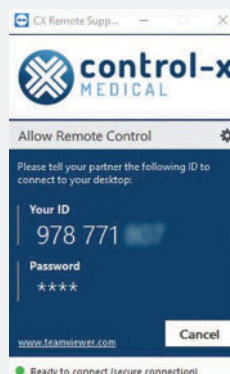

Please contact us at [sales@cxmed.com](mailto:sales@cxmed.com) or at +36 1 381 0301.

Zoomax Vetelev Gold Info sheet - Version 2.0 November 2021

## CERTIFICATIONS

Control-X Medical veterinary diagnostic products are CE marked.

